

Q-Plex™ Human Angiogenesis

The Q-Plex™ Human Angiogenesis (9-plex) is a fully quantitative ELISA-based chemiluminescent assay allowing the concurrent measurement of nine biomarkers or analytes (ANG-2, FGF basic, HGF, IL-8, PDGF-BB, TIMP-1, TIMP-2, TNFα, VEGF).

Each kit contains a 96-well plate, featuring the relevant biomarker panel, and all the reagents required to perform testing. Q-Plex plates are built by absorbing nine distinct capture antibodies in a defined array to the bottom of each well. Our high quality reagents help ensure the accuracy of your results.

Using just 50 µl of sample per well, up to 80 samples can be assayed for all nine markers in the panel within 2.5 hours. Q-Plex™ Arrays provide researchers with an easy-to-use and cost effective means of generating a cytokine profile for each sample.

Human Angiogenesis

Catalog # Product

150233HU Q-Plex™ Human Angiogenesis (9-plex)

Markers

ANG-2 FGF basic HGF IL-8
PDGF-BB TIMP-1 TIMP-2 TNFα
VEGF

Product Range

	Units	Calibrator Range*	LLD
Ang-2	pg/ml	10000 – 13.72	3.90
FGF basic	pg/ml	5000 – 6.9	12.7
HGF	pg/ml	14000 – 19.2	7.2
IL-8	pg/ml	2000 – 2.7	2.7
PDGF-BB	pg/ml	4000 – 5.5	2.2
TIMP-1	pg/ml	10000 – 13.7	10.2
TIMP-2	pg/ml	20000 – 27.4	15.9
TNFα	pg/ml	4000 – 5.5	5.2
VEGF	pg/ml	2000 – 2.7	2.7

* Actual values may vary from kit to kit. Please see the antigen card included in your kit for specific values.

Supplied by:



Tel: +44(0)1235 431390
sales@oxfordbiosystems.com
www.oxfordbiosystems.com

The Q-Plex™ Array

Spotting in 96 Well Plate

- ▶ Robotic liquid handlers print 20-50nL spots of capture antibody
- ▶ Each spot is a unique assay within the well
- ▶ Low spot-to-spot variability (CV)
- ▶ Spot size 350-500µm
- ▶ Plates are QC'd for spot quality

Performing the Assay

- ▶ Add as little as 50µL of sample
- ▶ Wash
- ▶ Add mix of detection antibodies specific to kit
- ▶ Wash
- ▶ Add streptavidin conjugated HRP or IR-Dye

Detection of Sample

- ▶ With the addition of substrate, a response is produced
- ▶ If antigen is present the spot emits a signal proportional to the amount of antigen in the sample
- ▶ If no antigen is present, the spot is not visible

Image Capture





- ▶ An image of the plate is taken via a high resolution camera contained within the Q-View Imager Pro or Q-View Imager LS
- ▶ The image file (TIFF) is imported into Q-View Software and a Q-View Project is created

Image Analysis

- ▶ Image is opened in Q-View Software
- ▶ Spots are automatically found on plate image
- ▶ Intensity of spot response is measured and raw data is generated
- ▶ User imports product specs and well layout

Data Analysis

- ▶ Raw data is analyzed and compared to in-plate standard
- ▶ Regression models used to calculate unknowns
- ▶ Standard curves are calculated and sample and statistical data is exported

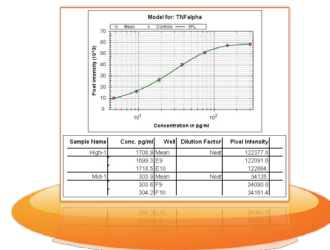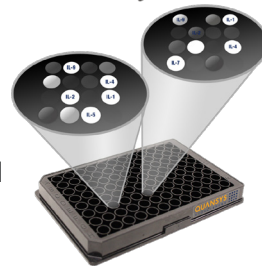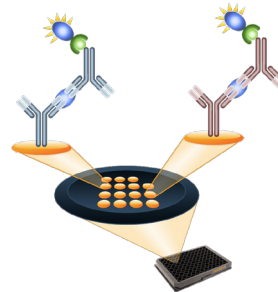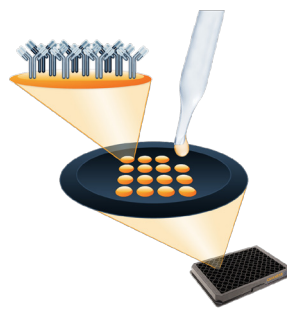

Increase productivity with Quansys products and services. Save Time, Sample and Money.

Q-Plex™ ARRAYS

Our Q-Plex Arrays are quantitative multiplex ELISAs with distinct proteins deposited in a defined array. Choose one of our standard kits for immediate delivery or customize the exact array you need.

Q-View™ IMAGERS

The Q-View Imager Pro and Q-View Imager LS are high quality, low-cost imaging systems for chemiluminescent assay imaging.

Q-View™ SOFTWARE

The Q-View Software is a user-friendly image analysis app that enables the acquisition and analysis of large amounts of multiplex ELISA data.

SAMPLE TESTING SERVICE

Our lab runs immunoassays on a wide variety of biological samples. Increase productivity and focus on your core research while we test your samples for you.



Quansys Biosciences
365 North 600 West
Logan, UT 84321
www.quansysbio.com

Fax: 435-750-6869
Phone: 435-752-0531
Toll Free: 888-QUANSYS (782-6797)

Quansys is an ISO 9001:2015 and ISO 13485:2016 registered company and complies with cGMP. Products are designed, developed, and manufactured according to the procedures outlined in our Quality Management System.

©2018 Quansys Biosciences Inc. All rights reserved. Q-View and Q-Plex are Trademarks of Quansys Biosciences Inc. ISO is a trademark of the International Standards Organization. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.

QUANSYS
B I O S C I E N C E S



www.quansysbio.com